


Installation instructions for ZEN 3.12 Smart Microscopy Workshop

- Extract *ZEN 3.12 installer.zip* files and run ZEISS microscopy installer via *setup.exe*
- Choose following options and install:

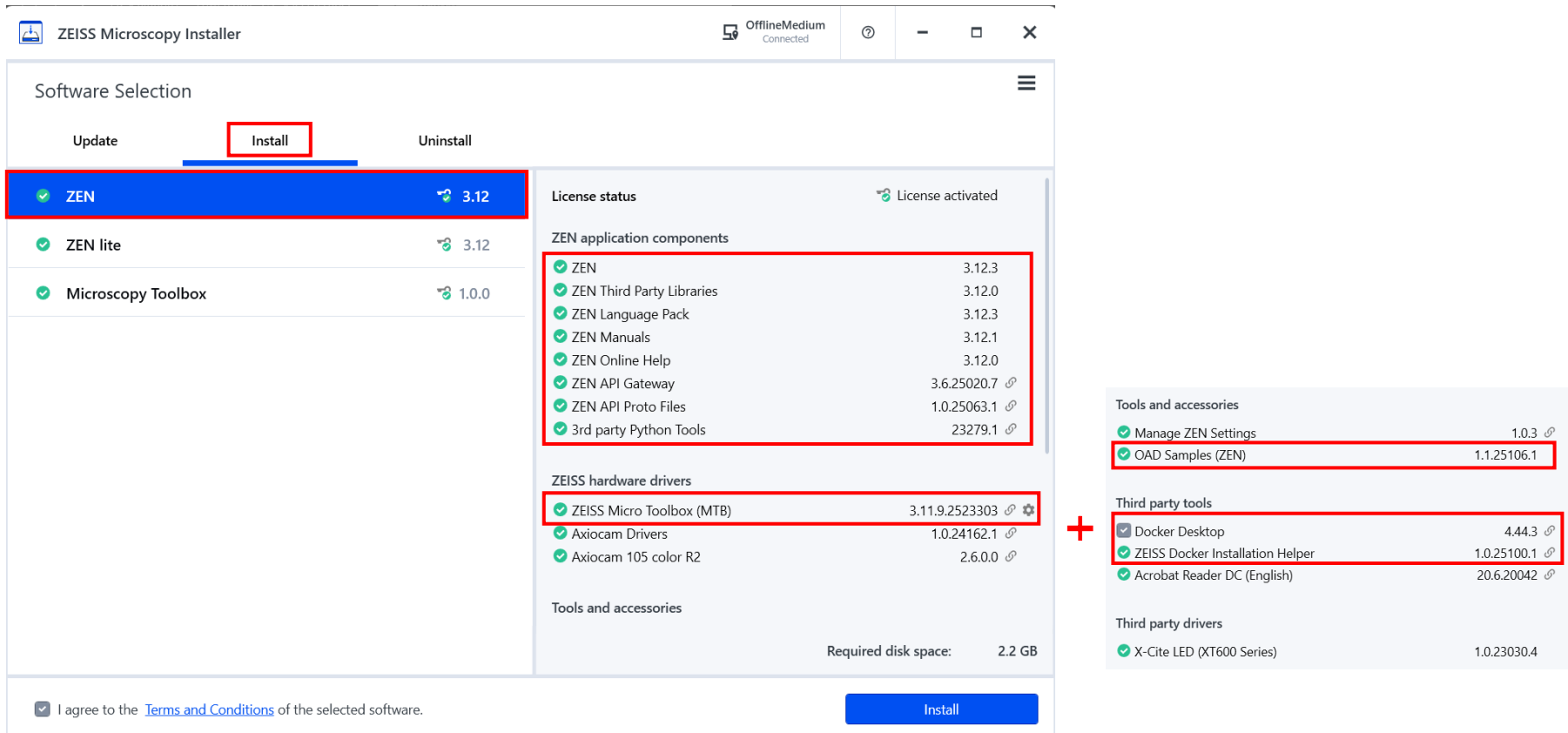

Software Selection

Update **Install** Uninstall

☒ ZEN 3.12

☒ ZEN lite 3.12

☒ Microscopy Toolbox 1.0.0

License status: License activated

ZEN application components

☒ ZEN	3.12.3
☒ ZEN Third Party Libraries	3.12.0
☒ ZEN Language Pack	3.12.3
☒ ZEN Manuals	3.12.1
☒ ZEN Online Help	3.12.0
☒ ZEN API Gateway	3.6.25020.7
☒ ZEN API Proto Files	1.0.25063.1
☒ 3rd party Python Tools	23279.1

ZEISS hardware drivers

☒ ZEISS Micro Toolbox (MTB)	3.11.9.2523303
☒ Axiocam Drivers	1.0.24162.1
☒ Axiocam 105 color R2	2.6.0.0

Tools and accessories

☒ Manage ZEN Settings	1.0.3
☒ OAD Samples (ZEN)	1.1.25106.1

Third party tools

☒ Docker Desktop	4.44.3
☒ ZEISS Docker Installation Helper	1.0.25100.1
☒ Acrobat Reader DC (English)	20.6.20042

Third party drivers

☒ X-Cite LED (XT600 Series)	1.0.23030.4
---	-------------

Required disk space: 2.2 GB

☒ I agree to the [Terms and Conditions](#) of the selected software.

Install